

# Freezerworks

## Basic Version Bar Code Package

The Freezerworks® Basic Bar Code Package allows you to print and read cryogenic bar code labels for your Freezerworks samples. **Here's what you get for \$2210:**



**Printer:** Zebra 2844® Thermal Transfer printer **specially built with a centered gap sensor - necessary for printing the Dataworks wrap-around labels.** Prints labels up to 4.09 inches wide at up to 4 inches per second. Self-centering label roll holder for simple label loading. 203 dpi print resolution. Will print all Freezerworks labels, using linear Code 128, DataMatrix 2D or PDF417. Smallest human readable print type is 23 cpi/6 points. Can be installed for a single workstation or as a networked printer.

**Keyspan® USB adapter:** necessary to print labels directly from Freezerworks 5.5



**Scanner:** Symbol DS6707® Digital Imager bar code scanner with attached USB cable. This scanner will read all codes printed through Freezerworks and Freezerworks Unlimited: DataMatrix 2D and PDF417 codes as well as Code 128 (Freezerworks standard and compressed). It is trigger and motion activated, and utilizes imager technology to read 2D barcodes from any rotation. It comes ready to read - just plug into a USB port.

**Intellistand:** Scanner stand - motion activated for hands-free scanning.



**Labels:** Dataworks wrap-around: excellent for cryogenic storage. The wrap-around design protects the label and maintains adhesion in ultra low temperature and multiple baths.


**1 ml vial size - 2.75" x 1.00" - fits most standard cryo vial (LBL-02WR)**

**OR**

**0.5 ml vial size - .625" x 2.125" - fits the smaller size vials (LBL-03WR)**

**Package includes 2400 labels and the necessary ribbon. Our standard is 1ml size unless specified otherwise.**

### Actual sizes shown

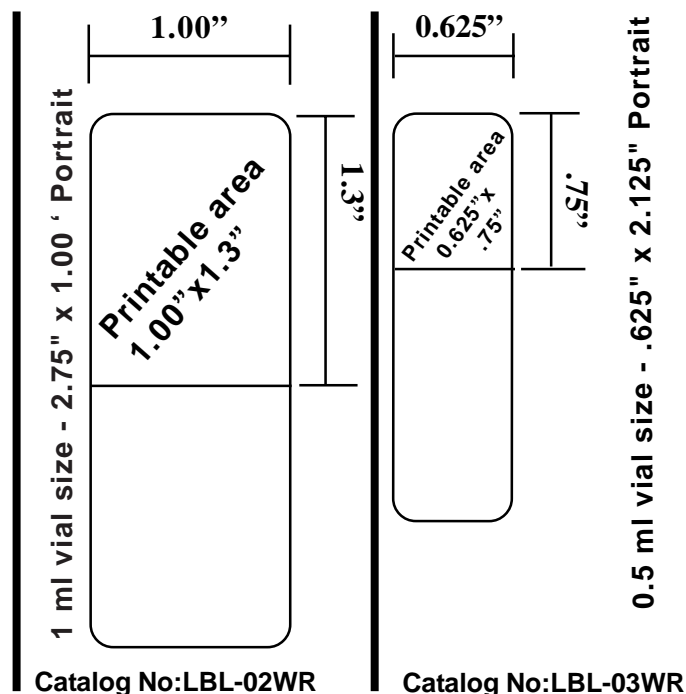

**Support:** Dataworks Customer Care Package: One full year of customer support and training via email, phone and web.